IMC

Morphecover

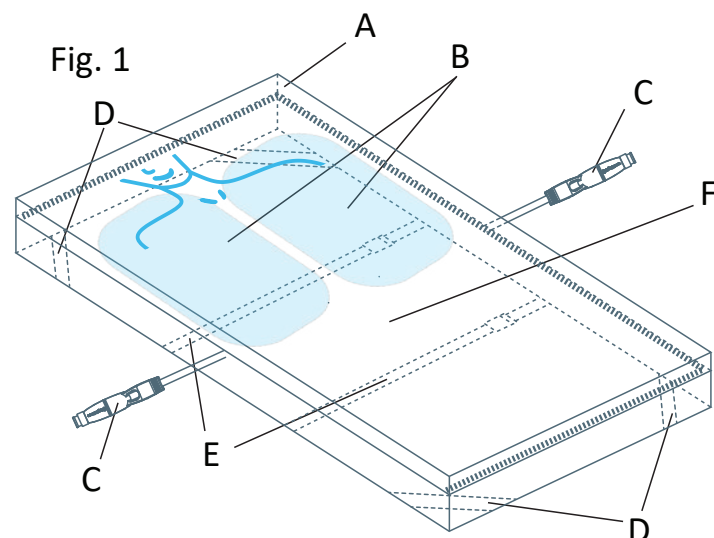
Fig. 2

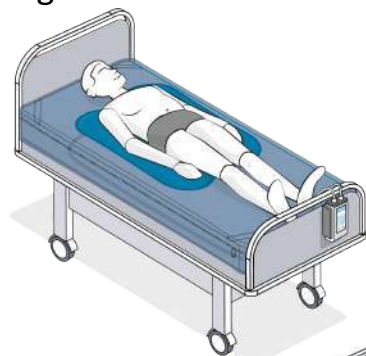
Fig. 3

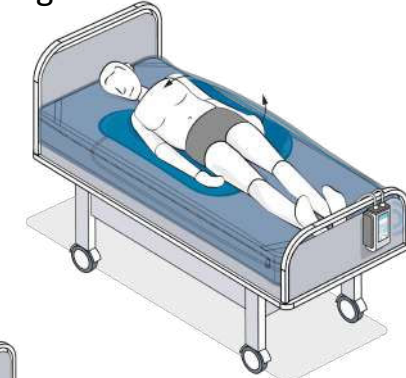
Fig. 4

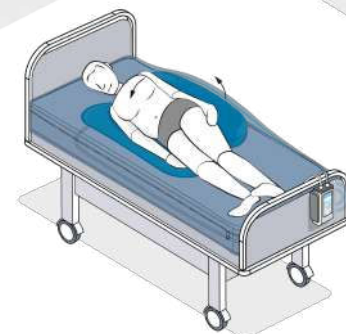
Fig. 5


5a

5b


5c

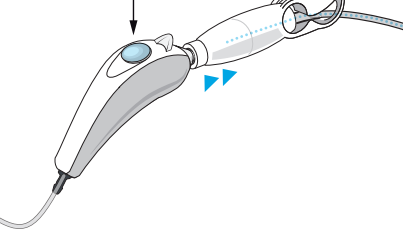5d


5e

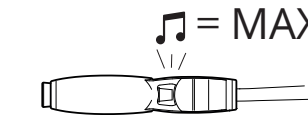5f

5g

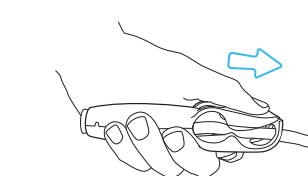Instructions for use (ENG) - Morphecover

Description

The Morphecover is a mattress cover with two inflatable cylinders. By inflating one of the cylinders using the Hummingbird (HMB) with compressed air, the patient is brought from supine in a stable lateral position (up to about 30° for a patient of 80kg/175lb and 20° for a patient of 180kg/400lb), both right and left.

The Morphecover consist of:

- Topper
- 2 inflatable cylinders
- 2 valves (detail in Fig. 5)
- 4 elastic corner straps
- 2 center straps
- Mattress cover

Intended purpose

Turning and maintaining a patient in a stable lateral position.

Installation instructions

- Place the Morphecover on the mattress with the silhouette to the head side of the bed.
- Place the elastic corner straps (D) around the underlying mattress.
- Lift your mattress and close the center straps (E) around the underlying mattress. Adjust the length of the center straps so that the Morphecover is well secured to the mattress.
- Apply bed linen.
- Check if the valves (C) are hanging on each side of the bed without bend or kink in the tube.

Operating instruction

- Position the patient on the Morphecover in supine position (Fig. 2).
- Check if the tube of the valve has no bend or kink and is not blocked by the bed. Inflate the Morphecover with compressed air using the HMB (Fig. 5). Hold the valve of the Morphecover with one hand while pushing the HMB against it with the other hand (Fig. 5c). The ring of the valve should be in closed position (Fig. 5a).

Note: When you hear the whistle at the start of the inflation, the tube is probably blocked. Recheck that the tube hangs freely and is not blocked by the bed.

- Stop inflating the Morphecover when the patient is in the desired position. Close the pinch valve (Fig. 5f).

Note: When you hear the whistle, the maximum filling has been reached.

- The patient will now remain in this lateral position. Check if the patient is laying correctly and doesn't touch the valves. Reposition the head and legs as per local procedure.

- To bring the patient back to supine position, push the pinch clamp forward with your thumb to unhook and open the lever (Fig. 5g). Slide the

- ring into open position to release the air from the Morphecover. (Fig. 5b.)
- Repeat step 2 to 4 if you want to tilt the patient to the other lateral side.

Warnings

The Morphecover is tested for 1000 inflation cycles. Leaks and deformation may occur after more usages.

Only inflate the Morphecover by using the Hummingbird (HMB) with compressed air at 4-6 bar/58-87psi. Higher pressures can cause damage. Stop inflating the Morphecover when the patient is at the desired position or when the overpressure valves open (whistle). Regularly check if patient is laying safe and comfortable.

Only for patients up to 180kg/400lb.

Make sure the bed side supports are high enough relative to the Morphecover.

Store under normal warehouse conditions.

Any serious incident related to the device must be reported to the manufacturer and the competent authority of the Member State where the user and/or the patient are/is established.

Cleaning instructions

Regularly/daily cleaning:

- First clean the product with water and soap. Use soft cloth and avoid harsh rubbing.
- Then, disinfect with isopropyl alcohol, chlorine dioxide or quaternary ammonium following the instruction of your disinfection product.

Disassembling instructions (for washing in the washing machine)

- Unlock and remove the valves (C) from the Morphecover.
- Open the zipper (located on the bottom side of the Morphecover) and remove the topper (A) from the Morphecover.
- Remove the cylinders (B) from the Morphecover.
- Keep the topper (A), cylinders (B) and valves (C) in a clean and dry place while washing the Morphecover.
- Close the zipper and the center straps (E) and clean the Morphecover in the washing machine (max 60°C/140°F). No drycleaning. Do not tumble dry.

Re-Installation after washing the Morphecover in the washing machine:

- Open the zipper and center straps (E) and place the cylinders (B) in the pockets. Place the tube of the cylinder through the opening of the pockets and trough the openings at the side of the Morphecover. The tubes should be at the bottom of the mattress cover (facing the ground).
- Place the topper (A) in the Morphecover and close the zipper.
- When the positioning of the topper in the Morphecover is not optimal, shake the Morphecover until the topper is correctly in place.
- Click the valves on the tubes of the cylinders.
- Follow steps of the installation guideline above.

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